User Manual

MULTIFUNCTIONAL

FACIAL SKIN CARE MASSAGER

Cleansing Rejuvenation Device

A facial beauty device with galvanic ion, heating, micro-vibration, and red light. This device aids in deep cleansing, rejuvenation, and facial lifting.

***Galvanic ionic:** Using the principle of a mild impulse current between the anode and cathode, it can soften and clean dirt and oil sebum deposits that lodge in the pores. This deep treatment acts as a natural cleanser that helps reduce pore size, and enriched charged oxygen ions activate nutrients to repair damaged skin.


***Micro Vibration:** Using physical vibration, and electric stimulation are both ways to increase elasticity and enhance cellular metabolism.

***Warm Massage:** Promotes absorption of skin-care products, helps makes skin texture tight and keeps the skin more elastic.

***Red Light:** 640nm-650nm red light deep into the dermis, helps stimulate new collagen and elastin growth and also reduces wrinkles and improves skin elasticity.


Cautions when using:

1. Press and hold the power button for 3 seconds to turn it on. It cannot be turned on when charging.
2. When using, hold the sensor strips on both sides of the wand, with the massage head in contact with the skin.
3. Please use it with a cotton pad when exporting the pores and remove the makeup remover or cleansing lotion with a cotton pad. Do not use it with facial cleanser/cleansing oil
4. This product uses micro-current ion technology. If you feel a small tingling sensation this is normal.
5. When using, do not stay in the same place for more than 3 seconds to avoid to skin irritation.

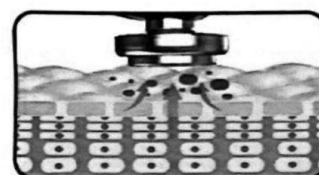
WARM CLEAN MODE

Warmth helps open the skin's pores, deep cleaning the skin to help electively refine the pores and aid in removing dark skin spots.

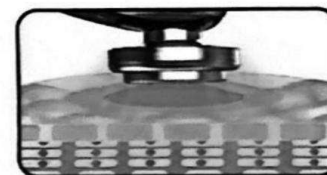PHOTOTHERAPY MODE

Drive skin-care products deep into the skin, activating the cell's elasticity. Using the red light helps aid in making the skin whiter and smooths out the wrinkles.

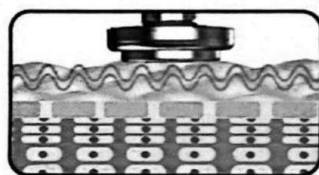SKIN LIFTING MODE

Promotes blood circulation, enhances skin elasticity and firming, and creates the best skin condition possible.

EYE CARE THERAPY MODE

Introduce eye cream or your favorite eye serum adequately, promote eye cream or eye serum absorption, achieve the maximum result of using your favorite eye care products.

Cautions:

1. Clean and rejuvenate, clean the dirt first and then use the cleanser (No need to use the device when cleaning face). Deep cleaning will pull the dirt from deep in the pores to the skin surface, and then wash the dirt with facial cleanser.
2. Each mode is timed for 4 minutes. After 4 minutes, it will enter standby mode. After 1 more minute, it will automatically shut down after no action.
3. When using, you need to hold the ion sensor strip on both sides of the fuselage and use the massage head to contact the skin to form a circuit. Please keep the skin moist when using otherwise, the device will not work properly.
4. If you use it on someone else's face, you need to contact the user's skin to form a circuit, otherwise, the device will not work properly.
5. Please use the device according to the length of time specified in the manual. Do not overuse it.

Product parameters

Product Name:	Photon Rejuvenation Wand
Model No:	W- 618
Product Size:	4.3*16.8cm
Product Weight:	About 90g
Package Contains:	1*device; 1*instruction manual, 1*USB cable
Input:	DC 3.7V/0.6A
Rated Voltage/ RatedCapacity:	3. 7V / 350mAH
Charging Time:	2-3 hours
Intensity:	3 levels
Country Of Origin:	China
Carried Standard:	GB4706.1-2005 GB-4706.15-2008

Note: In view of the consistent policy of continuous product improvement, the company reserves the right to improve product specifications, the actual product may vary.

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